

# *ProteinFind*<sup>®</sup> Goat Anti-Rabbit IgG(H+L), HRP Conjugate

Cat. No. HS101

Storage: at 2-8°C for one month; at -20°C for one year

## Description


Affinity purified *ProteinFind*<sup>®</sup> Goat Anti-Rabbit IgG(H+L) Antibody is a horseradish peroxidase (HRP) conjugated secondary antibody for ELISA and Western blot detection.

## Suggested Dilution

- Western Blot: 1:500-10000 dilution.
- ELISA: 1:500-5000 dilution.

## Notes

- Aliquot prior to use, and avoid repeated freezing and thawing.
- Diluted antibody should be used within 12 hours.
- Sodium azide is an inhibitor of HRP, thus avoid using sodium azide as a preservative in antibody dilution.

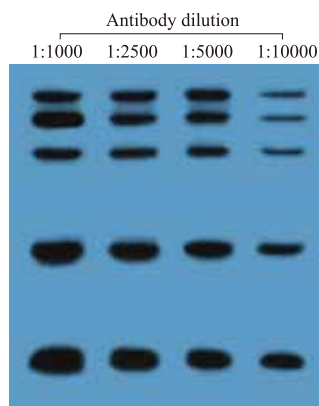

*EasySee*<sup>®</sup> Western Marker (5 µl/well)

**For research use only, not for clinical diagnosis.**

Service telephone +86-10-57815020

Service email [complaints@transgen.com](mailto:complaints@transgen.com)

