

ProteinFind[®] Goat Anti-Mouse IgG(H+L), HRP Conjugate

Cat. No. HS201

Storage: at 2-8°C for one month, at -20°C for one year

Description


Affinity purified *ProteinFind*[®] Goat Anti-Mouse IgG(H+L) Antibody is a horseradish peroxidase (HRP) conjugated secondary antibody for ELISA and Western blot detection.

Suggested Dilution

- Western Blot: 1:500-10000 dilution.
- ELISA: 1:500-5000 dilution.

Notes

- Aliquot prior to use, and avoid repeated freezing and thawing.
- Diluted antibody should be used within 12 hours.
- Sodium azide is an inhibitor of HRP, thus avoid using sodium azide as a preservative in antibody dilution.

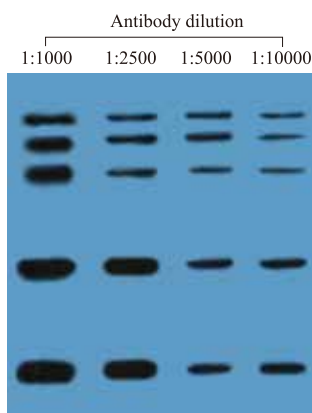

EasySee[®] Western Marker (5 µl/well)

For research use only, not for clinical diagnosis.

Service telephone +86-10-57815020

Service email complaints@transgen.com

