

Trans5K DNA Marker

Cat.No. BM141

Storage: at 2-8°C for six months; at -20°C for two years

Concentration: 0.095 mg/ml

Description

Trans5K DNA Marker is a premixed, ready-to-load molecular weight marker containing eight linear double-stranded DNA fragments. The DNA Ladder is suitable for use as molecular weight standards for agarose gel electrophoresis. The DNA Ladder contains 300 bp -5.0 kb DNA fragments. The 1.5 kb band (125 ng/5 μ l) has doubled intensity than other bands to serve as reference band.

Highlights

- Ready-to use molecular weight marker for DNA size range of 300 bp to 5.0 kb.
- All bands from restriction enzyme digested plasmids.
- Brighter 1.5 kb reference band.

Composition

300 bp, 500 bp, 800 bp, 1,000 bp, 1,500 bp (125 ng/5 μ l, the double intensity band), 2,000 bp, 3,000 bp, 5,000 bp


Storage Buffer

10 mM Tris-HCl (pH 8.4), 10 mM EDTA, 0.02% bromophenol blue, 5% glycerol

Notes

- In general, we recommend to use *Trans5K* DNA Marker per mm of gel lane, or 5 μ l per gel lane.
- The *Trans5K* DNA Marker is not intended for precise quantification of DNA mass but can be used for estimating the molecular weight of linear, double-stranded DNA fragments.
- Do not use the *Trans5K* DNA Marker with polyacrylamide gels.
- For long term storage, store at -20°C. If stored at -20°C, mix well after thawing.


Trans5K DNA Marker (5 μ l)

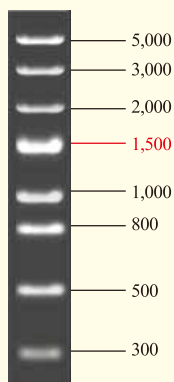

1.5% TAE agarose gel stained with GelStain
(loading volume: 5 μ l)

Recommend Agarose Gel Percentages for Trans DNA Markers

DNA Marker	Agarose	DNA Marker	Agarose
<i>Trans2K</i> [®] DNA Marker	1.5%	1Kb DNA Ladder	1.0%
<i>Trans2K</i> [®] Plus DNA Marker	1.5%	1Kb Plus DNA Ladder	1.0%
<i>Trans2K</i> [®] Plus II DNA Marker	1.5%	100bp DNA Ladder	2.0%
<i>Trans5K</i> DNA Marker	1.5%	100bp Plus DNA Ladder	2.0%
<i>Trans8K</i> DNA Marker	1.0%	100bp Plus II DNA Ladder	2.0%
<i>Trans15K</i> DNA Marker	0.8%	<i>Trans</i> DNA Marker I	2.0%
		<i>Trans</i> DNA Marker II	2.0%

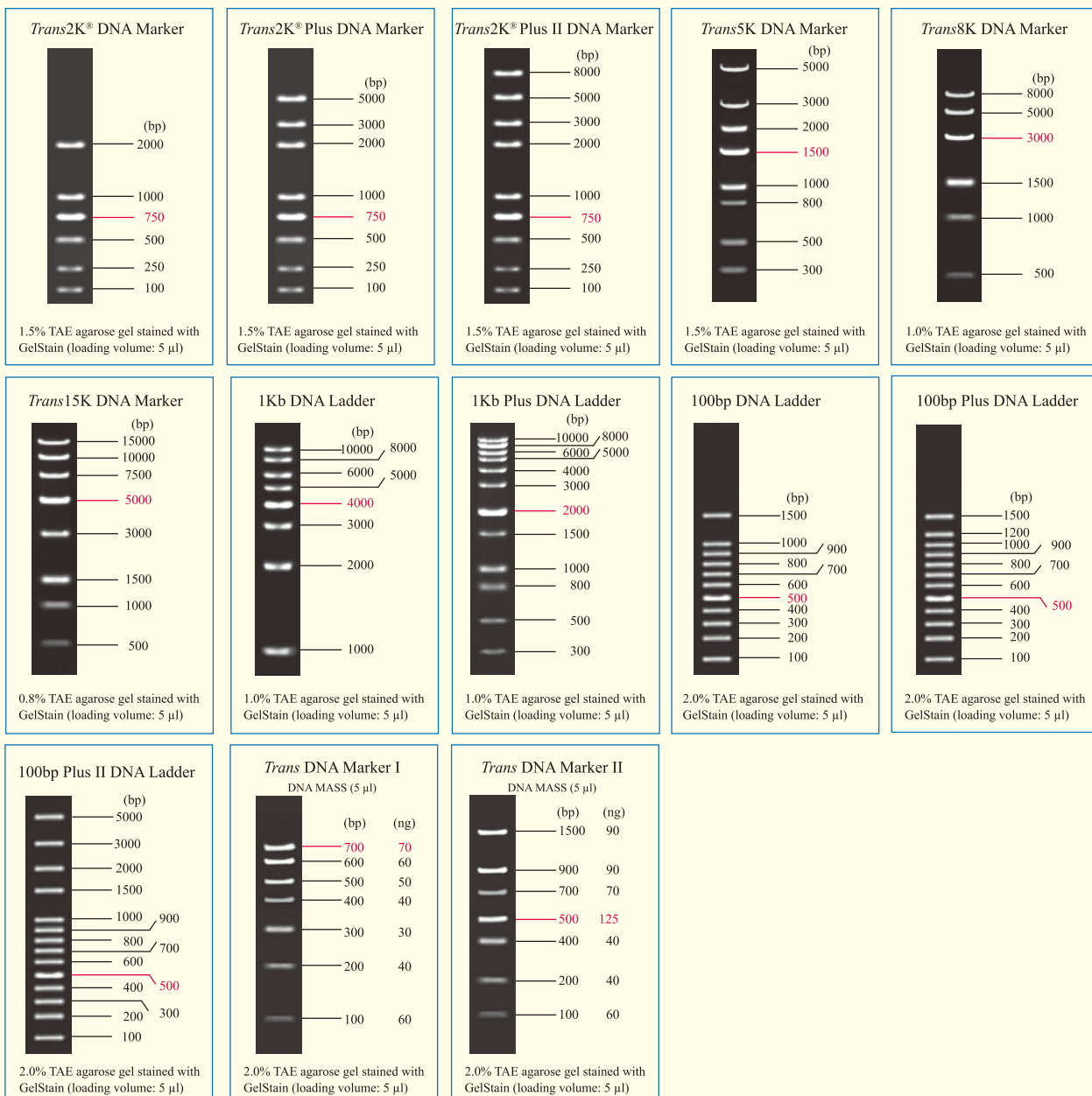

FOR RESEARCH USE ONLY